



BEERE MEDICAL® TRU-CONTROL® RATCHETING DRIVERS

AT A GLANCE

- Ergonomically-designed driver has one of the world's strongest, most precise ratcheting mechanisms.
- Tru-Control® ratcheting mechanism is tested to 200 in-lbs of torque.
- Smooth, precise ratcheting with minimal play.
- A variety of adapters are available. See chart on opposite side.

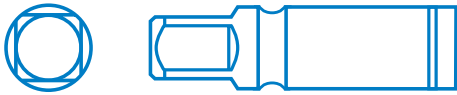
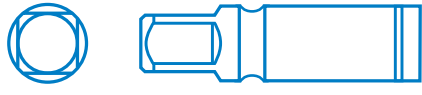

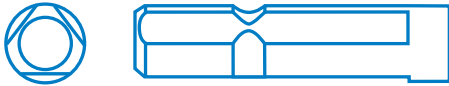




PRECISE AND STRONG TRU-CONTROL® RATCHETING MECHANISMS

- Tru-Control® ratcheting driver is tested to 200 in-lbs of torque.
- Proprietary design provides smooth, precise ratcheting with minimal radial, axial, or lateral play.
- Low back-torque makes it easier to start/insert screws.
- You can ratchet in forward or reverse direction, or lock in a central position for standard driving applications.
- Mechanism is designed to provide a smooth, consistent feel.
- Ratcheting mechanism can be housed in a variety of Ortho-Grip® handles and mated with various standard connectors.
- Outstanding reliability: Ratcheting mechanism has a lifetime guarantee on the materials and workmanship.*

ERGONOMIC ORTHO-GRIP® HANDLES

- Teleflex Medical OEM offers a large selection of handle platforms that are compatible with most surgical instruments.
- Ergonomic Ortho-Grip® silicone handles are:
 - Autoclaveable
 - Crack-resistant
 - Texturized for a non-slip grip
 - Designed to fit existing trays
- Custom configure Ortho-Grip® handles from many design options, which include:
 - Cannulated or solid designs
 - Impaction caps
 - External QuickLock® adapters
- Choose virtually any Pantone® color for the handle.**

ADAPTER OPTIONS AND SHAFT CONFIGURATIONS

Adapter Options	Male Shaft Configurations	Design Standards***
1/4" Square QuickLock® Quick-Lock® adapters prevent the accidental release of the shaft from the handle.		DES-002 (DES-0061: spring-loaded ball)
3/16" Square QuickLock®		DES-0011 (DES-0013: spring-loaded ball)
A-O		DES-0701
Tri-Flat or Tri-Lobe		DES-0800
Hex (6mm)		DES-0903
Hudson, Small		DES-1101
Hudson, Large		DES-1102
Jacob-Hall (Zimmer-Hall)		DES-1201

***Design standards are available to assist with your handle configuration project.

**ORTHO-GRIP®
HANDLE PLATFORMS**

The Tru-Control® ratcheting mechanism is available on many Ortho-Grip® handle platforms, including:

- Axial: Small, mini, and large
- Contour: One, two, or three scallops
- T-Handle: Standard and comfort
- Palm
- Pin Puller
- Scallop

TELEFLEX MEDICAL OEM

5307 95th Avenue
Kenosha, WI USA

Toll Free: 800.295.8505 (USA only)

Phone: 262.657.2800

Fax: 262.657.2801

E-mail:

oeminfo@teleflexmedical.com

Website:

www.teleflexmedialoem.com

Imagine it. Achieve it.®

©2010 Teleflex Incorporated. All rights reserved. Teleflex, Beere Medical, Tru-Control, Ortho-Grip, QuickLock, and "Imagine It. Achieve it." are registered trademarks of Teleflex Incorporated or its affiliates. 512-1010

*Guarantee is void if instrument was abused, used in a manner other than its intended purpose, or improperly maintained.

**Pantone is a registered trademark of Pantone, LLC.